

Instructions

How to Use




- 1. Short press the power button to select the desired mode.
- 2. Press and hold the power button for more than 0.5 seconds to start operation after following step 1.
- 3. If you want to change modes during operation, short press the power button to switch.



- 4. Place the nozzle in your mouth and use your lips to prevent water from splashing around.

※ Use in front of the sink to allow used water to fall into the sink.

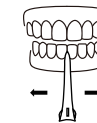

- 5. Position the nozzle tip horizontally.
- 6. Clean by moving along the gumline from one end to the other.

Instructions




- 7. Move the nozzle to clean both the inner and outer sides of the teeth.




- 8. If you want to turn off the power during operation, press and hold the power button for more than 0.5 seconds. It will automatically turn off after two minutes of use or when the water is depleted.

※ If the hose is tightly attached to the water tank wall, the power may not turn off, even if the water is depleted.

Detachment Method

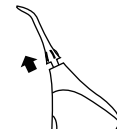

- 1. Pull the nozzle vertically to detach it. Pulling in the wrong direction may cause damage.

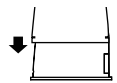

- 2. When detaching the water tank, pull it vertically.

※ Since there is no separate locking mechanism, you can detach it directly

Precautions

Precautions

- Do not expose to flames or heat.
- Do not disassemble or modify.
- Batteries with a visibly damaged appearance may be dangerous, so contact the manufacturer or seller immediately.
- Specifications may change without notice for the purpose of performance improvement.

During Use

- Familiarize yourself with the instructions and precautions before use.
- If you have recently had gum surgery or experience pain, consult with a professional before use.
- Initial use may cause gum bleeding due to the pressure. If bleeding occurs, it is recommended that you use the soft mode. If bleeding persists for more than two weeks, consult a professional.
- If your gums swell significantly after use, discontinue use and consult a professional.
- Do not use on eyes, ears, nose, or other body parts.
- Keep away from children.
- Continuous operation for more than five minutes may cause motor overheating and malfunction.
- If used upside down, water may leak through the air hole in the water tank. This is not a malfunction, and you can continue to use it.
- Avoid strong impacts or excessive force.
- Do not put foreign substances or hot water into the water tank. Doing so may cause a malfunction or harm to the user.

During Charging

- Do not place the charger in a damp place.
- Dry the charging terminal thoroughly before charging.
- Use only the designated charger.
- Do not allow metal objects to touch the electronic terminals.
- Do not charge the battery at temperatures above 40 degrees Celsius.
- If storing the battery for a long time after charging, remove the charger after it is fully charged.
- Charge the battery only in the specified manner.
- The charging time may vary depending on the battery level, temperature, and environment.

Cleaning and Maintenance

- When replacing the nozzle, use the dedicated nozzle provided or sold separately. The device is not compatible with nozzles from other manufacturers.
- After using the product, clean the nozzle and let it air dry. If it is not completely dry, bacteria may breed.
- Once three months have passed since use or if the nozzle is deformed, replace the nozzle.
- If visible damage to the magnet is confirmed, replace the nozzle.
- Do not immerse in boiling water or chemicals for disinfection, as doing so may cause a malfunction.
- Do not wash in devices such as dishwashers or washing machines for cleaning.

Caution

FAQ

Symptoms	Causes	Solutions
Product not working	Battery discharge	Please charge before use.
Leakage between nozzle and main body	Installation error	Check for foreign substances between the nozzle and the attachment, then reattach.
Leakage from the water tank	Installation error	Check for foreign substances on the water tank and its connection, then reattach.
Installation error	Reverse usage	Water may leak through the air holes if used upside down. This is not a malfunction, so feel free to continue using it.
Decreased output	Low battery	Please charge before use.
	Nozzle deformation	Replace the nozzle if there is any foreign substance buildup or deformation on the nozzle or attachment.
	Incorrect mode setting	The mode might be set to soft. Change the mode before use.
No water spray during operation	Mesh filter clogged	Clean the mesh filter.
	Empty water tank	Fill the water tank.
Unable to charge	Faulty charging connector	If there is any foreign substance such as water or dust on the charging connector terminal, clean it thoroughly before use.
	Product malfunction	If charging doesn't work even after cleaning the charging port, request after-sales service.

Battery Information

·This product supports Type-C charging.

·During charging, the battery indicator will be red if the battery level is 80% or less and blue if it's more. During use, it will be red if the battery level is 20% or less and blue if it's more.

Battery Level	Battery Indicator
Charging, over 80%	Blue
Charging, below 80%	Red blinking
In use, over 20%	Blue
In use, below 20%	Red

· The charging specification for this product is 5 V 1 A. Please use an appropriate adapter. (Charging may not work with adapters below 5 V 1 A.)

· Ensure there is no moisture when connecting USB power. Moisture may cause fire or electric shock.

Product Specifications

Product Information

- Product Name: SONIQ WATER FLOSSER
- Model: SWF-001~010
- Material: ABS, PC, Silicone, etc.
- Rated Input: DC 5 V
- Power Consumption: 5 W
- Charging Time: 5 hours
- Pulse Pressure: 1,800/minute
- Battery Specification: 3.7 V, 2,550 mAh
- Waterproof Rating: IPX7
- Size: 91 x 73 x 227 (mm)
- Weight: 395 g

After-Sales Service and Maintenance

The free after-sales service period is six months. Even within the free period, repairs due to consumer-caused faults will be subject to charges. For inquiries about after-sales service and maintenance, refer to the manufacturer's information.

Manufacturer

- Country of Manufacture: South Korea
- Manufacturer: Genoss Co., Ltd.
- Address: D-Factory, 56 Changryongdaero 256beon-gil, Yeongtong-gu, Suwon-si, Gyeonggi-do, Republic of Korea
- Tel. 031-888-5100
- Website: www.genoss.com

Disposal

- Do not dispose of the battery recklessly.
- Dispose of the product according to the disposal method regulations of each country and region.




SONIQ
WATER FLOSSER


Model: SWF

Dentium GENOSS


Components

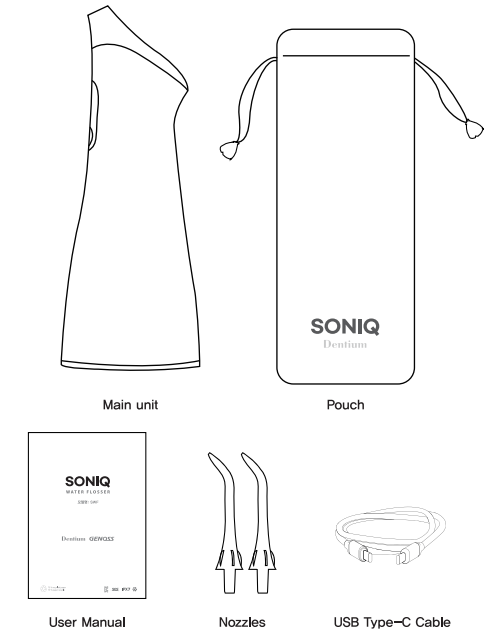


Product Structure

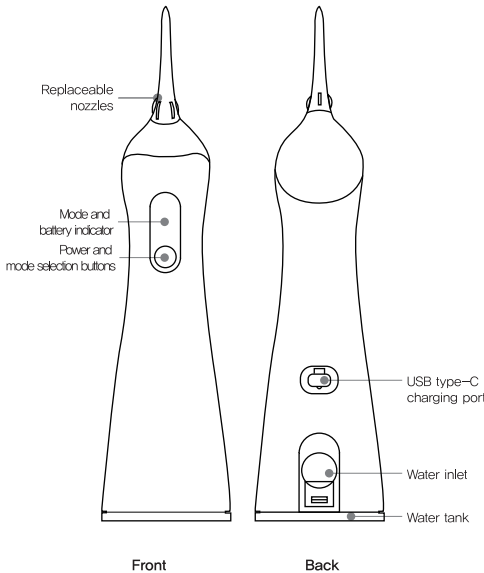

Product Features

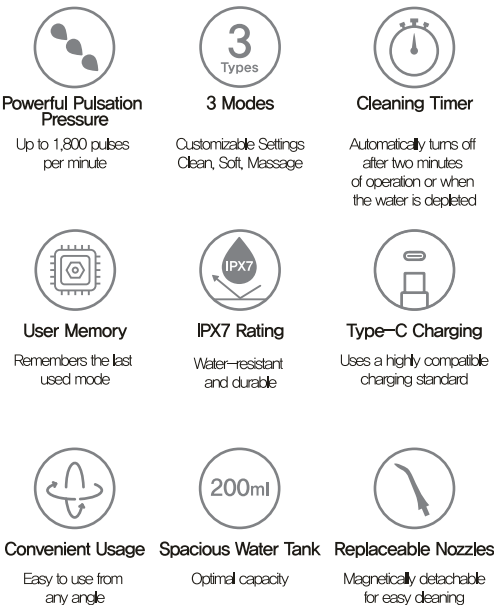

Usage Mode Explanation

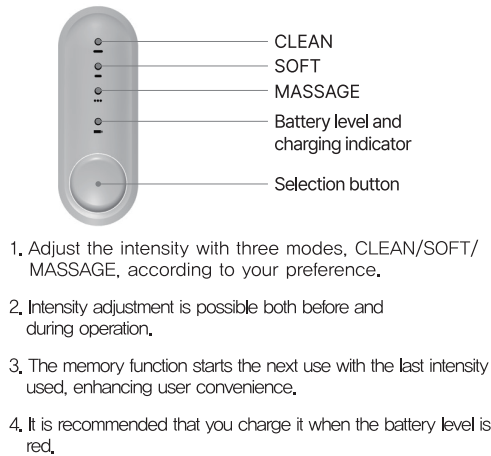


Usage Modes


Mode	Mode Description
CLEAN	Strong water pressure mode
SOFT	Gentle water pressure mode
MASSAGE	Mode that automatically adjusts the intensity for a gentle gum massage

Usage

Charging Method

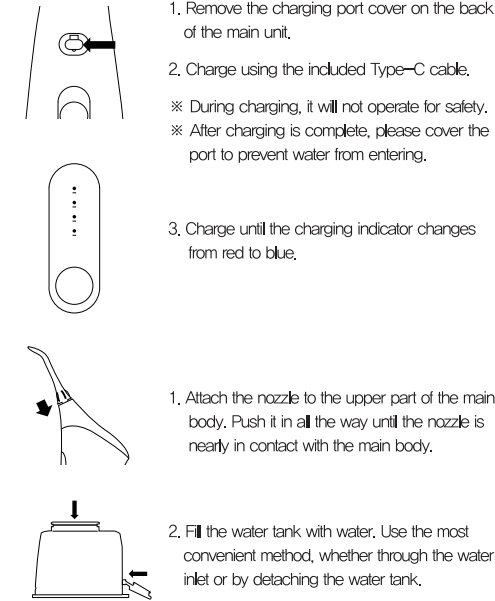

Storage and Cleaning

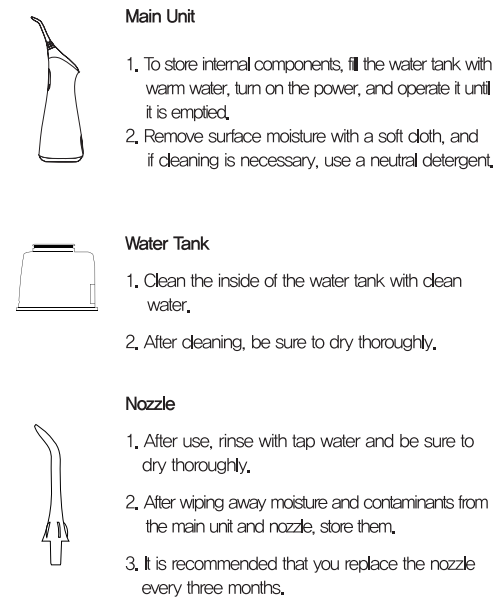

Caution Do not expose any components to low temperatures with moisture after use.